




50bp DNA Ladder

For research use only.

Catalog No.	Specification	Storage/Shelflife
EDL106-01	250 μ l	-20C°/2 years
EDL106-02	500 μ l	-20C°/2 years

Introduction

The 50bp Ladder DNA Marker is composed of 9 DNA fragments with a length of 50 bp to 500 bp, dissolved in 1 \times Loading Buffer, and can be used for direct electrophoresis of 5-10 μ l during used, which is very convenient to use.

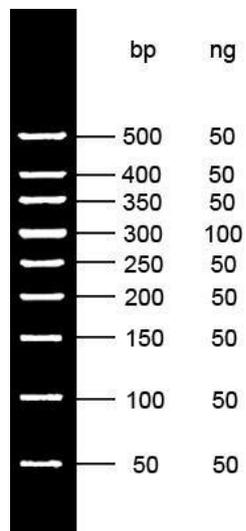

Note

1. When the width of the sample well during electrophoresis is less than 5 mm, a clear band can be obtained by taking 5 μ l Marker electrophoresis each time. If the sample hole is widened, the sampler volume of Marker must be increased appropriately.
2. For DNA electrophoresis, the purity of Agarose greatly affects the clarity of DNA bands. Therefore, the best quality Agarose should be used for electrophoresis. The recommended gel concentration is 2.5% to 3%.
3. When performing Agarose electrophoresis, the concentration of Agarose is closely related to the separation performance of DNA fragments. The greater the concentration of Agarose, the better the separation performance of short DNA; conversely, the smaller the concentration of Agarose, the more favorable for the separation of long DNA.

Gentaur Molecular Products BVBA

Address: Voortstraat 49, 1910 Kampenhout, Belgium

T: 0032 16 58 90 45 | E: info@gentaur.com

Websites: www.gentaur.com | www.maxanim.com