

### MAD345Hu22

Monoclonal Antibody to Casein Kinase 1 Delta (CSNK1d)

Organism Species: Homo sapiens (Human)

Instruction manual

FOR RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)

# [PROPERTIES]

**Source:** Monoclonal antibody preparation

Host: Mouse

Antibody isotype: IgG2b Kappa

**Purification:** Protein A + Protein G affinity chromatography

Clone number: C1

Traits: Liquid

Concentration: 1mg/ml

**UOM:** 200µl

**Cross Reactivity: Porcine** 

Applications: WB; IHC; ICC; IP.

### [ IMMUNOGEN ]

Immunogen: Recombinant CSNK1d (Met1~Lys409) expressed in E.coli

Accession No.: RPD345Hu01

#### [ APPLICATIONS ]

Western blotting: 0.5-2µg/mL;

Immunohistochemistry: 5-20µg/mL;

Immunocytochemistry: 5-20µg/mL;

Optimal working dilutions must be determined by end user.

#### [FORMULATION]

Form & Buffer: Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300,

50% glycerol.

### [ STORAGE AND STABILITY ]

Storage: Avoid repeated freeze/thaw cycles.


Store at 4°C for frequent use.

Aliquot and store at -20°C for 24 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined

by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

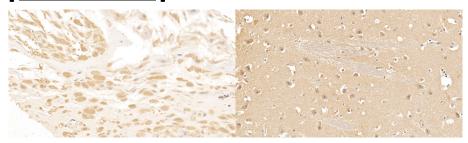

Human Placenta Tissue; Primary Ab: 10μg/ml Mouse Anti-Human CSNK1d Antibody Second Ab: 2μg/mL HRP-

DAB staining on IHC-P; Samples:

Linked Caprine Anti-Mouse IgG Polyclonal Antibody (Catalog:

SAA544Mu19)

DAB staining on IHC-P;

Samples: Human Cerebrum Tissue;

Primary Ab: 10µg/ml Mouse Anti-

Human CSNK1d Antibody

Second Ab: 2µg/mL HRP-Linked

Caprine Anti-Mouse IgG Polyclonal

Antibody

(Catalog: SAA544Mu19)

### [ IMPORTANT NOTE ]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.

Address: Voortstraat 49, 1910 Kampenhout, Belgium

T: 0032 16 58 90 45 | E: info@gentaur.com

Websites: www.gentaur.com | www.maxanim.com